

# GreatSmile Breeze

## Enjoy your new Great Smile Breeze – Oral Water Jet.

(Please read this instruction leaflet carefully).


**EASY** – with the flick of the switch and 30sec washing adds greatly to your oral health

**CONVENIENT** – only 30sec while in the shower and eliminates splash and mess that you get from counter top oral irrigator systems

**EFFECTIVE** – better results than brushing and flossing alone

### How to install your Great Smile Breeze:

- 1 Unscrew showerhead (A) from shower pipe/hose (B).
- 2 Wrap plumber's tape (included) once or twice clockwise around shower pipe/hose threads (B). Pull tight enough to draw tape down into threads.

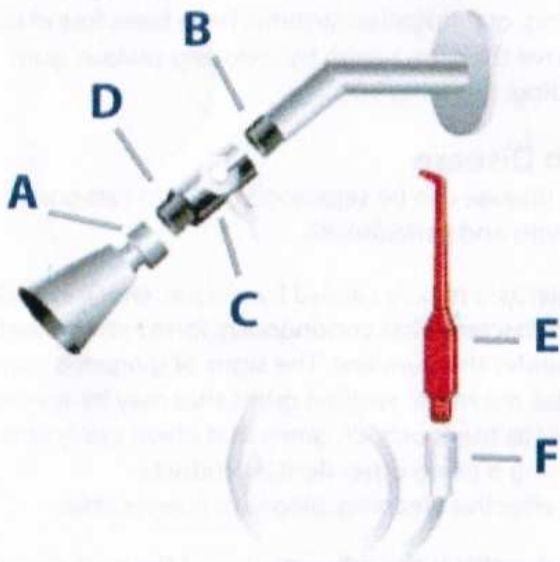

- 3 Screw *Great Smile Breeze* (C) onto shower pipe/hose clockwise by hand until tight. Do not use wrench.
- 4 Wrap plumber's tape around *Great Smile Breeze* threads (D) and screw showerhead onto threads and hand tighten.
- 5 Attach hand-held nozzle (E) to connector at end of *Great Smile Breeze* tube (F) by pressing and turning until it locks into place. Do not force.
- 6 Turn on shower and use selection switch to choose *Great Smile Breeze* or shower and to adjust water pressure.
- 7 Twist and pull nozzle (E) out of hose connector (F) and store in holder (on wall).

**TIP:** To eliminate tube coiling, run hot water through the hose and hold the hose straight and tight for several seconds.

### How to use your Great Smile Breeze effectively:

**Switch On** - Set selection switch to *Great Smile Breeze* and breeze away.

**Breeze** - Ensure you hold your *Great Smile Breeze* close to the gums, almost touching them.

- Wash your teeth and gums for 30 seconds, pausing between each tooth as you slowly move along the gum line.

- Continue until all areas between and around your teeth have been thoroughly irrigated.

**Switch Off** - Set selection switch back to shower

*Great Smile Breeze* should be used with a low water pressure system, and works most effectively using warm water only. Prior to using your *Great Smile Breeze* in your mouth, test the water pressure, by holding it toward the shower wall. By starting with a low water pressure, you can increase or adjust the water pressure to your personal preference.

Initially, it is normal to feel a slight tickling sensation, which will soon diminish as your gums become accustomed to regular stimulation. If you should experience pain or unusual discomfort, lower the pressure setting right away. Always follow the specialist recommendations of your dental professional.

